

DONJOY[®]

Clima-Flex[™] OA

More Flexibility,
Better Mobility

NOW AVAILABLE TO PRESCRIBE IN THE UK

 **DJO**[®]



DONJOY®

CLIMA-FLEX™ OA KNEE BRACE

More Flexibility, Better Mobility

Clima-Flex™ OA offers a new level of comfortable, cooling relief through compression and offloading for patients with mild to moderate OA knee pain. New DonJoy C-6 Fabric Technology delivers breathable, moisture-wicking comfort that's soft to the touch, antimicrobial and incorporates temperature control.

Keyless offloading allows quick and easy hinge adjustment, offering patients convenient, comforting relief from OA knee pain. Clima-Flex™ OA is ideal for everyday or light leisure activities such as walking, golf, or hiking.

- **Clima** = Temperature Control
- **Flex** = Mobility & Comfort

KEYLESS, DUAL DIAL OFF LOADING

New, easy-to-use integrated dual-action dial allows quick and easy hinge adjustment for comfortable, convenient offloading of the knee anytime, anywhere.

Keyless Dual Dial Telescoping hinge provides extra 3° of off-load, in addition to the 3° built into the frame.

E-Z WRAPAROUND APPLICATION

E-Z Grip handles and intuitive wraparound design ensure form-fitting comfort. Brace compression level is easily adjusted based on patient needs and comfort.

PERFECT FIT & CONTROL

Single upright with semi-rigid thigh and calf cuffs lined with foam, provide comfortable support.






NEW

C-6 Technology™ Advanced Cooling Material

Pure carbon threads are woven into this moisture-wicking, highly breathable and cooling antimicrobial fabric to naturally repel odour and provide temperature control.

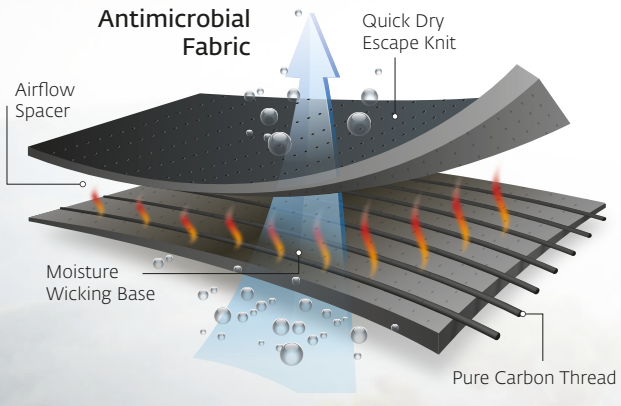

WHY BRACE FOR KNEE OSTEOARTHRITIS?

Unloading Knee Braces are Cost-Effective for Those Awaiting Surgery

A 2017 UK study suggested that unloading knee braces are a cost-effective method as a bridging therapy and long-term use can delay and reduce surgical intervention. It also provided evidence that if patients wear the brace for more than 2 years, surgical intervention is not needed.¹

Quoted in the latest NICE Guidelines for Osteoarthritis Care and Management

“People with osteoarthritis who have biomechanical joint pain or instability should be considered for assessment for bracing/joint supports/insoles as an adjunct to their core treatments.”²

Size	Thigh Circumference 15cm Above Mid-Patella	Knee Centre Circumference Mid-Patella	Calf Circumference 15cm Below Mid-Patella
1 (XS)	33 - 39 cm	30.5 - 33 cm	25.5 - 30.5 cm
2 (S)	39 - 47 cm	33 - 35.5 cm	30.5 - 35.5 cm
3 (M)	47 - 53 cm	35.5 - 38 cm	35.5 - 40.5 cm
4 (L)	53 - 60 cm	38 - 43 cm	40.5 - 47 cm
5 (XL)	60 - 67 cm	43 - 48.25 cm	57 - 50.75 cm
6 (XXL)	67 - 75 cm	48.25 - 53.25 cm	50.75 - 56 cm
7 (XXXL)	75 - 81 cm	53.25 - 58.5 cm	56 - 61 cm

Size	Left Leg Medial Compartment or Right Leg Lateral Compartment	Right Leg Medial Compartment or Left Leg Lateral Compartment
	PIP Code	PIP Code
1 (XS)	4185690	4185757
2 (S)	4185658	4185773
3 (M)	4185682	4185765
4 (L)	4185708	4185740
5 (XL)	4185666	4185732
6 (XXL)	4185674	4185716
7 (XXXL)	4185641	4185724

Measurement and Sizing

- Measure at three points on the leg and use the chart to determine the size required.
 - Thigh 15cm above mid-patella (kneecap)
 - Knee mid-patella (kneecap)
 - Calf 15cm below mid-patella (kneecap)
- If one of the three measurements is a size up, select that size. For example, Thigh Size 2, Knee Size 3, Calf Size 2. Order Size 3

Note that the off-loading hinge is located on the opposite side of the affected compartment.

References

1. Lee PYF, Winfield TG, Harris SRS, et al. Unloading knee brace is a cost-effective method to bridge and delay surgery in unicompartmental knee arthritis. *BMJ Open Sport Exerc Med* 2017;2:e000195.doi:10.1136/bmjsem-2016-000195
2. www.nice.org.uk/guidance/cg177



Watch the fitting video